PURPOSE / OVERVIEW

This guide provides instructions for installing the ASI Floating Tray Assembly Kit for both the Glider (90-2148) and the GliderLite (90-1148). Follow these steps during initial setup of your ASI dental delivery system.

PARTS INCLUDED

Glider (90-2148):

Description	Qty
Floating Tray	1
8" Arm Assembly	2
12" Arm Assembly	2
Riser Assembly	2

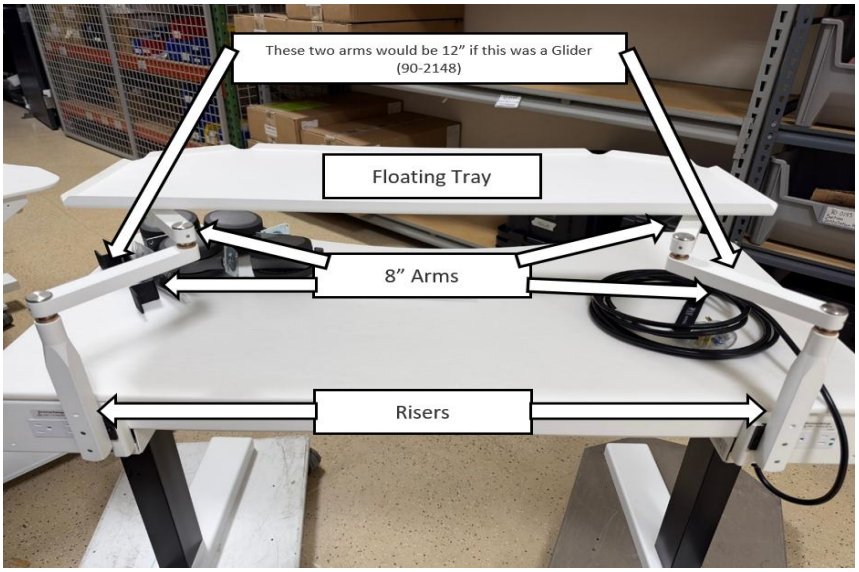
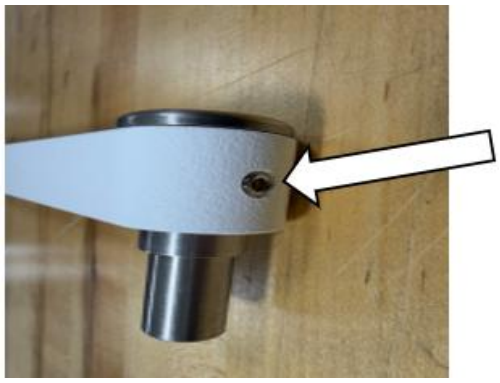
GliderLite (90-1148):

Description	Qty
Floating Tray	1
8" Arm Assembly	4
Riser Assembly	2

TOOLS REQUIRED

- 3/32” Hex Allen Key (as needed)

INSTRUCTIONS

1	• Slide the risers into the dovetail bracket on the back of the cart. • Insert the arm shafts into the tops of the risers. • Connect a second arm to each installed arm. • Attach the arms to the floating tray. ***NOTE**: The Glider (90-2148) will have the risers, then 12" Arms, and then the 8" Arms, and then the tray.*  **Fig. 1: GliderLite Floating Tray**
2	The arms can be tightened as needed with a 3/32" Allen Hex Key in the set screws on the back of the arms.  **Fig. 2**

QUESTIONS OR ASSISTANCE

For technical support, contact ASI Dental Specialties:

Web: **ASIsupport.com** | Email: **ASIservice@ASIsupport.com** | Phone: **303.407.6072** | Parts: **ASI.parts**

REVISION HISTORY (For Internal Master Part File Only)

CO #	Rev	Description of Change	Author	Approval Date
CO-26-187	A	Initial Release	P. Wick	10/21/26